

極メカ

[KIWAMIMECHA]

Extreme mechanisms

Exquisite comfort that approaches the deep, powerful and pliant touch of real hands

With improved shape of the knead balls and drive gear pitch, the extreme mechanisms allow extending the maximum thrust distance to 100 mm. High performance 32-bit CPU reproduces delicate and versatile knead techniques. The positioning function allows adjusting the knead ball position at 11 steps back and forth and vertically, so as to massage at exact points of stiffness.

AIR BAGS

Soft and powerful massaging with the air bags with 35 air chambers

A large capacity and high output air compressor quickly drives the air bags arranged all over the body. Each air bag is designed to fit a particular body section.

Fuji's original SOLUTION MESSAGE that combines the "MECHA. & AIR" MESSAGES.

Solution Massage is Fuji's synthetic healing technology that offers a total massage for relaxation and comfort by integrally controlling the "MECHA. & AIR" MESSAGE.

Depending on body sections and degrees of stiffness or strain, gentle surface stimulation by air and strong point stimulation by mechanical movements are applied in combination.



MECHA. MESSAGE

Severely stiffened sections are massaged intensively with stimulation by the knead balls.

- NECK
- SHOULDER
- SPINE
- WAIST

AIR MESSAGE

Soft surface stimulation by the air bags massage in a manner to push muscles up.

- SHOULDER
- ARM
- BACK
- WAIST
- HIP
- THIGH
- CALF
- LEG
- SOLE

7-STEP ADJUSTMENT

Strength or gentle is adjustable depending on the mood or body condition at the moment!

5-STEP ADJUSTMENT

※3-step adjustment for the arm

MECHA. × AIR

A total 559 types of combination massage

More versatile and delicate massaging is achieved by combining the MECHA. & AIR-MESSAGES with the basic patterns of EXTREME MECHA., MECHA. & AIR-MESSAGE.

EXTREME: 5 types
MECHA.: Basic 50 types
AIR: Basic 9 types
MECHA. × AIR: 495 types


Total: 559 types

※ All illustrations represent images just for explanation.

ORIGINAL FUNCTION

Motor-driven sliding foot-rest

Motor-driven foot-rest can slide back and forth in a range of approx. 180 mm. Adjust it at the optimum position when using the chair.

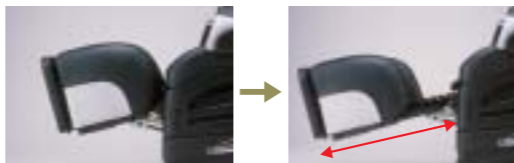

Maximum adjustable length, approx. 180mm

MOTOR-DRIVEN RECLINING

Reclining angle can be adjusted at desired positions from the remote controller at hand.



[KIWAMIWAZA]
New-極技



BK (Black)



CS (Beige)

CYBER-RELAX

Massage Chair EC-3000

■Measurement (approx.): (With foot-rest stored) 76W × 123D × 122H (cm) [29.9 × 48.4 × 48.0 inches]. (with foot-rest reclined) 76W × 201D × 80H (cm) [29.9 × 79.1 × 31.5 inches]

■Weight (approx.): 86kg (190LB) ■Reclining angle: 120 - 170°

■Rated time: 30 min. ■Power consumption: 115W (50/60 Hz)



REMOVABLE BACK AND SEAT PADS, AND PILLOW

Each part can be removed easily for maintenance. ※Photograph shows the black type.



SOLE PAD (ACCESSORY)

Antibacterial specification sole pads for use on the bottom of the foot-rest.



Optional item



Side table AS-TS2



SAFETY PRECAUTIONS

- Be sure read the "User Manual" carefully beforehand to use the product properly.
- The following persons are prohibited to use the product, as it could cause physical problems.
 - Persons who are prohibited to receive massage, for example: those suffering from thrombosis (embolism), severe aneurysm, acute varicose veins, various kinds of dermatitis and skin infections (including inflammations of subcutaneous tissue), etc.
 - The following persons must consult a doctor before using the product, as its use could cause physical problems:
 - Those who suffer from malignant tumors.
 - Those who are in an unstable condition in early pregnancy or immediately after childbirth.
 - Those who suffer from anesthesia related to peripheral circulation trouble due to diabetes, etc.
 - Those who have a body temperature higher than 38°C (pyrogenic period) (Example: When strong acute inflammatory symptoms are observed, such as the malaise, chill, anisopiasis, etc. When the person feels weak.)
 - Those with wounds in the locations to be treated.
 - Those who suffer from acute illness (dolorific)
 - Those who suffer from osteoporosis.
 - Those whose spines are broken.
 - Those who suffer from inflammatory diseases such as sprains, muscle strain, etc.
 - Those who suffer from a sensory disorder.
 - Those who have a heart disease.
 - Those who use a medical electronic device embedded in the body such as a pace-maker, which is susceptible to electromagnetic interference.
 - Those who have some abnormality in the spine or whose back is bent to the right or left.
 - Those who suffer from anemia.
 - Those who need to massage a part of body treated medically in the past.
 - Who need repose.
 - Those who are under a doctor's treatment.
 - Those who feel any abnormality on the body other than the above.

Fuji Medical Instruments Mfg.Co., Ltd.

● For further information;

* Design and specifications are subject to change without prior notice for improvement of the product.
* Colors of printing and actual product may differ slightly.

FUJIIRYOKI

Relieves mental and physical strain gradually, with a sense like being enwrapped.



New-極技 [KIWAMIWAZA]

CYBER-RELAX

Massage Chair EC-3000

Delicate fingertip touch and pliant movements of hands and arms that detect exactly the massage points, which vary by person. We have developed professional techniques that remove accumulated strain, even in your brain!

